



Abcb4-KO Genotyping Protocol

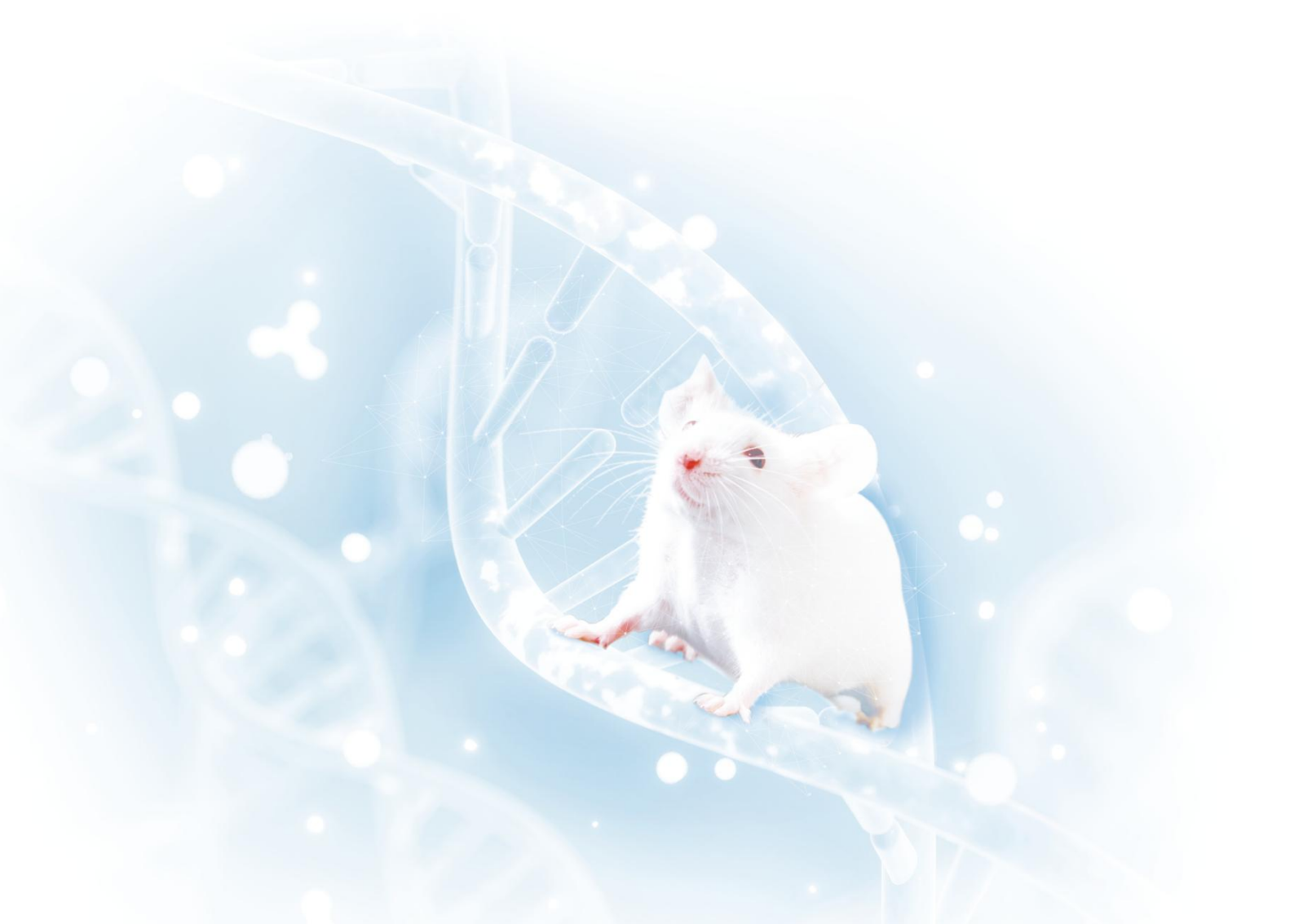

Common Name	Abcb4-KO	Cat. NO.	NM-KO-241600
Strain of Origin	C57BL/6J	Version	V1

Genotyping strategy

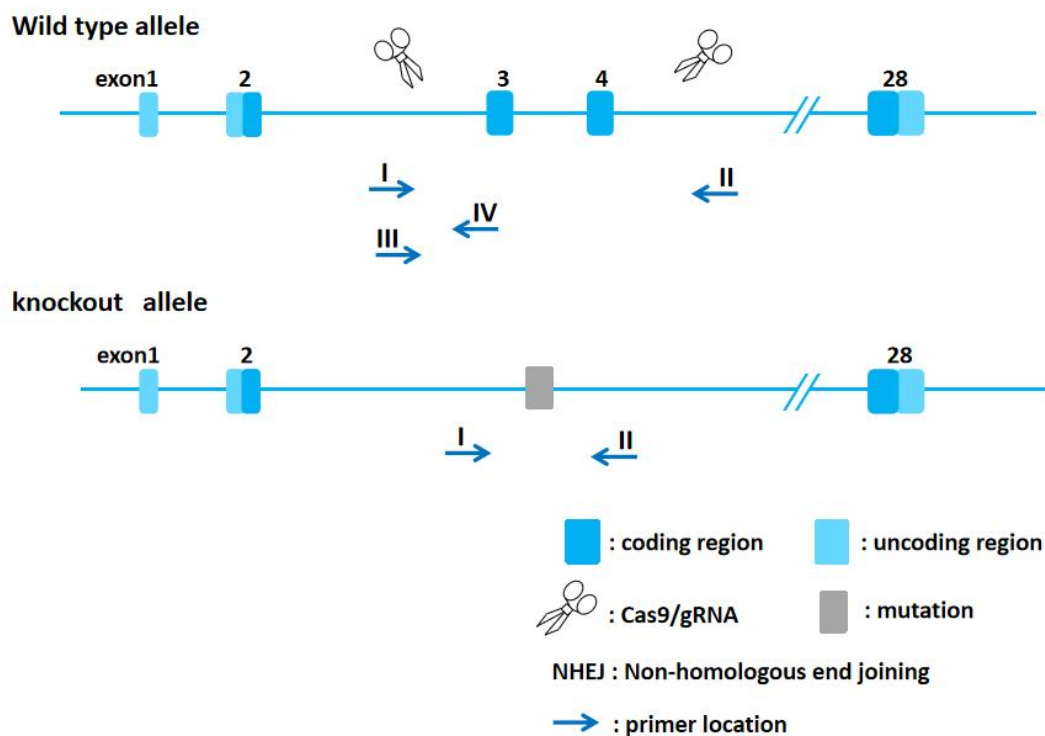

Primers

Primer	Sequence (5'→3')	Primer type
P1	CCCTCAGAGAGTGTGGACTG	Forward
P2	AACCCAGACAGAAGCCTACAA	Reverse
P3	TCCAGAAAAGCAGTCTCCCTAT	Forward
P4	TACCA GTGTCAACAGGCCAA	Reverse

Expected results

Results	
Genotype	Knockout type: -3465bp Wild type: P1P2 =4121 bp; P3P4 =396 bp Heterozygote: P1P2 =4121 bp and 656 bp; P3P4 =396 bp Homozygote: P1P2 =656 bp

Note: In both wild-type and heterozygous mice, whether the P1 and P2 primers can amplify larger bands does not affect the interpretation of the results, because the purpose of designing this pair of primers is to amplify KO band

Reaction &Cycling

PCR Reaction System	Reaction Component			Volume (μl)
	ddH2O			8.0
	2×Rapid Taq Master Mix			10.0
	P1(10 pmol/μl) or P3(10 pmol/μl)			0.5
	P2(10 pmol/μl) or P4(10 pmol/μl)			0.5
	Genomic DNA			1.0
	Total			20
	2×Rapid Taq Master Mix from Vazyme(Code Number: P222-01)			
Cycling Reaction	Step	Temp	Time	Note
	1	95°C	5 min	
	2	95°C	20 sec	
	3	60°C	20 sec	
	4	72°C	20 sec	35 repeats to 2
	5	72°C	5 min	
	6	12°C	Hold	